


Toileting Sling & Toileting Sling High



Padded Arms
Reinforced Back
Divided Legs

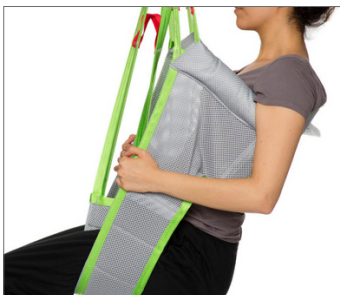

Our Toileting Slings feature a wide, divided leg partition that enables fast and simple application for dressing and undressing during toileting or bathing. They are suitable for patients with poor trunk stability and feature padding under the arms for extra comfort and are available with extra head support or without.

Three Options

These Toileting slings are available in two low back models and one high back model with adjustable neck support for increased head stability. The wide partition of the leg parts makes the sling application fast and simple and enables dressing and undressing during toileting or bathing. All models can be used with or without an optional waist belt that is available in an assortment of sizes with a Velcro or buckle closure.

Optimal Positioning

The low back models feature an upright seated position that makes transferring to and from a toilet much easier. The high back model features a reclined (backward tilted)

position that prevents the user from sliding out of the sling and is especially beneficial for those hypertonic users that have a tendency to overstress themselves (spastic performance) or those hypotonic users with lesser body control.

Easy Care and Easy Fit

Our slings are manufactured to the highest quality standards. They feature reinforced straps and high strength nylon webbing around the sling to make them strong enough to resist even the most active use. Each sling label is colour coded to ensure quick and easy selection of the correct sling and the straps are colour coded for easy determination of required positioning.



Walking Aids



Healthcare Beds



Lifting Solutions



Hygiene



Convertible chairs

Accessories



Extension Loops
Used to extend lift straps in order to alter lifting position or to hold leg parts together. Assorted lengths. Part no. 21180.*



Waist Belts
Velcro or buckle enclosure. Available in Junior and Adult sizing. Velcro Part no. 25145-** Buckle Part no. 25146-***



Upholstered Cover
Protects user from hanger bar impact. Not for use in a wet environment. Part no. 50866



Hanger Bars for 2-, 3- or 4-point Suspension
Part no. 57551, 55764 [2- or 3-Point - Quick Release], 55619 [4-Point - Attach to leg hook]

Lifting Situations



Use to/from a bed
Suitable for transfers and repositioning to and from a bed. Can be left behind patient for short period of time.



Use to/from a chair
Suitable for transfers and repositioning to and from a chair. Can be left behind patient for short period of time.




Use for showering
Only Polyester Net fabric can be used for this application.




Use for toileting
Only Polyester Net fabric can be used for this application.


Care and Features

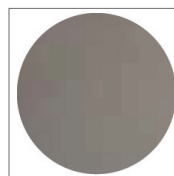

Polyester
Durable, low-friction, easy to clean, apply and remove.

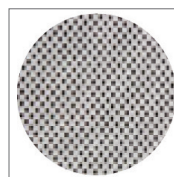

Polyest Net
Soft, thin material that breathes and dries quickly. Ideal for bathing or showering.

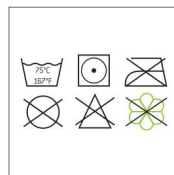

Care
Care instructions located on sling label. Use a hospital grade laundry bag while laundering.



Technical Specifications

Part Number			Size	Colour	Model #25080/25085			Model #25090			Model #25080/25085	Model #25090
Toileting Sling High Net	Toileting Sling Net	Toileting Sling Polyester			A	B	C	A	B	C		
-	25085-JS	-	Junior Small	JS	28cm / 11"	62cm / 24.4"	22cm / 8.7"	-	-	-		
25090-JM	25085-JM	25080-JM	Junior Medium	JM	28cm / 11"	67cm / 26.3"	26cm / 10.2"	28cm / 11"	73cm / 28.7"	70cm / 27.6"		
25090-S	25085-S	25080-S	Small	S	32cm / 12.6"	74cm / 29.1"	30cm / 11.8"	32cm / 12.6"	82cm / 32.3"	74cm / 29.1"		
25090-M	25085-M	25080-M	Medium	M	48cm / 18.9"	90cm / 35.4"	33cm / 13"	48cm / 18.9"	92cm / 36.2"	83cm / 32.7"		
25090-L	25085-L	25080-L	Large	L	49cm / 19.3"	97cm / 38.2"	33cm / 13"	49cm / 19.3"	99cm / 39"	90cm / 35.4"		
25090-XL	25085-XL	25080-XL	Extra Large	XL	56cm / 22"	100cm / 39.4"	38cm / 15"	56cm / 22"	101cm / 39.8"	100cm / 39.4"		
Maximum User Weight			300 kg / 660 lb									

* Available in assorted lengths: 10cm/3.9" (21180-10), 15cm/5.9" (21180-15), 20cm/7.8" (21180-20), 25cm/9.8" (21180-25), 35cm/13.7" (21180-35)

** Junior Sizes: 40cm/15.7" (25145-40B), 50cm/19.7" (25145-50B), 60cm/23.6" (25145-60B); Adult Sizes: 60cm/23.6"/(25145-60), 70cm/27.6" (25145-70), 80cm/31.5" (25145-80), 90cm/35.4" (25145-90), 110cm/43.3" (25145-110), 120cm/47.3" (25145-120)

*** Junior Sizes: 40cm/15.7" (25146-40B), 50cm/19.7" (25146-50B), 60cm/23.6" (25146-60B); Adult Sizes: 60cm/23.6"/(25146-60), 70cm/27.6" (25146-70), 80cm/31.5" (25146-80), 90cm/35.4" (25146-90), 100cm/39.4" (25146-100), 110cm/43.3" (25146-110), 120cm/47.3" (25146-120)



We improve freedom of movement and quality of care

www.humancaregroup.com

HUMAN CARE SWEDEN (HQ)
Årstaängsvägen 21B
117 43 Stockholm
+46 8 665 35 00
info@humancaregroup.com

HUMAN CARE CANADA
10-155 Colonnade Road
ON K2E 7K1 Ottawa
+1 613 723 6734
info.ca@humancaregroup.com

HUMAN CARE THE NETHERLANDS
Elsperweg 124
8076 PA Vierhouten
+31 577 41 21 71
info.nl@humancaregroup.com

HUMAN CARE UNITED STATES
8006 Cameron rd, Suite K
TX 78754 Austin
+1 5 124 76 7199
info.us@humancaregroup.com